

## INTRODUCTION

The INCOGNITO™ T7 Ψ-RNA Transcription Kit\* is optimized for high-yield synthesis of pseudouridine-containing RNA (Ψ-RNA) from an *in vitro* transcription (IVT) reaction. Although yield varies with the DNA template and other factors, a standard 2 hour, 20 μl reaction will yield up to 100 μg of RNA from 1 μg of the control template. These yields are made possible by the high-performance properties of the T7-Scribe enzyme. The standard reaction can be scaled up to produce milligram amounts of RNA containing the canonical nucleotides ATP, CTP, GTP and the modified nucleotide pseudouridine-5'-triphosphate (ΨTP).

It has been shown that Ψ-mRNAs are translated into protein at higher levels and induce lower innate immune responses in human and other mammalian cells that express various RNA sensors compared to corresponding unmodified mRNAs.<sup>1-4</sup>

INCOGNITO T7 IVT Ψ-RNA can be further processed into mRNA (5'-capped and 3'-poly[A] tailed). Options include: post-transcriptional capping using CELLSCRIPT's ScriptCap™ Cap 1 Capping System (contains both ScriptCap Capping Enzyme and 2'-O-Methyltransferase); post-transcriptional tailing using CELLSCRIPT's A-Plus™ Poly(A) Polymerase Tailing Kit; and/or co-transcriptional tailing through the use of a template-encoded tail.

## MATERIALS


### Materials Supplied

**Important** Store at –20°C in a freezer without a defrost cycle. Do not store at –70°C.

INCOGNITO™ T7 Ψ-RNA Transcription Kit Contents (10 reactions)	
Kit Component	Reagent Volume
T7-Scribe™ Enzyme Solution	20 μl
10X T7-Scribe™ Transcription Buffer	20 μl
GAΨC PreMix (25 mM each GTP, ATP, ΨTP & CTP)	60 μl
100 mM Dithiothreitol (DTT)	20 μl
RNase-Free DNase I, 1 U/μl	10 μl
ScriptGuard™ RNase Inhibitor, 40 U/μl	10 μl
T7 Control Template DNA, 0.5 μg/μl	10 μl
RNase-Free Water	1.4 ml



\* see patent and license information on page 9.

 For more information, consult the appropriate safety data sheet (SDS) at [www.cellscript.com/products.html](http://www.cellscript.com/products.html).

**T7 Control Template DNA:** Is a linearized 4.1-kb plasmid that contains a T7 promoter followed by a phage lambda dsDNA insert that encodes a 1,375-base runoff transcript. The Control Template DNA is provided at a concentration of 0.5 µg/µl in T<sub>10</sub>E<sub>1</sub> Buffer (10 mM Tris-HCl, pH 7.5, 1 mM EDTA).

### Materials Required, but not Supplied

- A DNA template for transcription of your RNA of interest
- Materials or kits for purification of the RNA product  
(For suggestions, see Section C "Purification of the Ψ-RNA Transcription Product")
- RNase-free TE Buffer (10 mM Tris-HCl, pH 7.5, 1 mM EDTA)

## SPECIFICATIONS

### Storage Buffers

RNase-Free DNase I is provided in a 50% glycerol solution containing 50 mM Tris-HCl, pH 7.5, 10 mM CaCl<sub>2</sub>, 10 mM MgCl<sub>2</sub> and 0.1% Triton<sup>®</sup> X-100. ScriptGuard RNase Inhibitor is supplied in a 50% glycerol solution containing 50 mM Tris-HCl, pH 7.5, 100 mM NaCl, 10 mM DTT, 0.1 mM EDTA and 0.1% Triton X-100. All other enzymes are provided in a 50% glycerol solution containing 50 mM Tris-HCl, pH 7.5, 100 mM NaCl, 1 mM DTT, 0.1 mM EDTA and 0.1% Triton X-100.

### Unit Definition

One unit of RNase-Free DNase I digests 1 µg of pUC19 DNA to oligodeoxynucleotides in 10 minutes at 37°C.

One unit of ScriptGuard RNase Inhibitor results in 50% inhibition of 5 ng of RNase A. Activity is measured by the inhibition of hydrolysis of cyclic 2',3'-CMP by RNase A.

### Functional Testing

The INCOGNITO T7 Ψ-RNA Transcription Kit is functionally tested under standard reaction conditions using the T7 Control Template DNA. The kit must produce at least 90 µg of RNA from 1 µg of the T7 Control Template DNA in 2 hours at 37°C.

### Contaminating Activity Assays

All components of the INCOGNITO T7 Ψ-RNA Transcription Kit are free of detectable RNase and DNase activity, except for the inherent activity of the RNase-Free DNase I component.

**BEFORE YOU START: IMPORTANT TIPS FOR OPTIMAL *IN VITRO* TRANSCRIPTION****◆ Template Requirements:**

The optimal templates for *in vitro* transcription are linear double-stranded DNA (dsDNA) molecules with 5'-protruding ends. DNA templates with blunt ends are less preferable and **templates with 3'-protruding ends should not be used.**

Transcription templates can be prepared from clones of the DNA to be transcribed in plasmids or other circular dsDNA vectors by linearizing the vectors downstream of the cloned DNA using a suitable restriction endonuclease or other means.

Alternatively, transcription templates can be generated by PCR amplification of RNA or DNA of interest using a strategy that results in joining of a T7 promoter to the appropriate end of the PCR product (e.g., wherein, the T7 promoter is either joined to the DNA or RNA that is amplified or is incorporated into one of the PCR primers).

**◆ Template Efficiency and Incubation Time:**

*In vitro* transcription of 1 µg of the T7 Control Template DNA using the INCOGNITO T7 Ψ-RNA Transcription Kit yields approximately 100 µg of ~1.4-kb RNA in 2 hours at 37°C in a standard 20 µl reaction.

However, yields vary for different templates based on the template sequence, structure, length, purity and the sequence and length of the particular RNA polymerase promoter. Examples of contaminants that can affect transcription yield include RNase, phenol, trace metals and SDS. See the Technical Appendix for suggestions related to template purification.

**◆ Optimizing the Reaction:**

The recommended reaction conditions should give excellent yields of RNA with most templates.

However, the yield may be improved for some templates by extending the reaction time (e.g., to 4 hrs), increasing the amount of template in the reaction, or increasing the reaction temperature from 37°C to 42°C.

**◆ Maintaining an RNase-Free Environment:**

Highly stable RNases are ubiquitous, including on human skin.

Creating an RNase-free work environment and maintaining RNase-free solutions is critical for successful work with RNA.

We strongly recommend to:

- Use RNase-free tubes and pipette tips.
- Always wear gloves when handling kit components or samples containing RNA and change gloves frequently, especially after touching potential sources of RNase contamination such as door knobs, pens, pencils and human skin. Do not touch any kit component or tube containing RNA with an ungloved hand.
- Keep all kit components tightly sealed when not in use. Keep all tubes containing RNA tightly sealed during the incubation steps.

## PROCEDURE

## A. Synthesis of IVT RNA

1. Set up the transcription reaction **at room temperature** by adding the reagents **in the order indicated below**:

Standard INCOGNITO T7 Ψ-RNA IVT Reaction	
Component	Amount
RNase-Free Water	x μl
Linearized template DNA with T7 RNAP promoter	1 μg
10X T7-Scribe Transcription Buffer	2 μl
GAΨC PreMix	6 μl
100 mM DTT	2 μl
ScriptGuard RNase Inhibitor	0.5 μl
T7-Scribe Enzyme Solution	2 μl
Total Reaction Volume	20 μl

**Important** Assemble transcription reactions at room temperature in the order indicated at left. Assembly of transcription reactions at <22°C or in an alternate order, can result in the formation of an insoluble precipitate.



Transcription Buffer stored at -70°C may result in the formation of a white precipitate. To dissolve it, heat the tube at 37°C for 5 minutes and mix thoroughly.



One microgram of linear dsDNA template is recommended for most reactions. If the concentration of the DNA template is <0.2 μg/μl, concentrate it and then resuspend in the appropriate amount of RNase-Free Water.

2. Incubate at 37°C for 2 hours.

## B. DNase I Treatment of IVT Reaction

1. DNase I treatment is used to remove the DNA template from the IVT reaction.

Standard DNase I Treatment of Transcription Reaction	
Component	Amount
IVT Reaction (from Step A)	20 μl
RNase-Free DNase I	1 μl
Total Reaction Volume	21 μl

2. Incubate for 15 minutes at 37°C.
3. Proceed to RNA Purification.

### C. Purification of the Ψ-RNA Transcription Product

Purify the RNA using your preferred method. The method chosen should remove residual proteins and unincorporated NTPs from the RNA. Several options are listed below. RNA can be stored at  $-20^{\circ}\text{C}$  or  $-70^{\circ}\text{C}$ . For long-term storage, RNA can be stored as an ethanol pellet.

I) **Ammonium Acetate Precipitation:** Selectively precipitates RNA, while leaving most of the protein and unincorporated NTPs in the supernatant. Note: for this method, the RNA to be purified must be >100 bases in size.

- 1) Add one volume of 5 M ammonium acetate (21  $\mu\text{l}$  for the standard reaction), mix well.
- 2) Incubate for 15 minutes on ice.
- 3) Pellet the RNA by centrifugation at  $>10,000 \times g$  for 15 minutes at  $4^{\circ}\text{C}$ .
- 4) Remove the supernatant with a pipette and gently rinse the pellet with 70% ethanol.
- 5) Remove the 70% ethanol with a pipette without disturbing the RNA pellet.
- 6) Allow pellet to dry, then resuspend in RNase-Free Water, TE or other suitable buffer.
- 7) While usually unnecessary, steps 1-6 may be repeated a second time for even cleaner RNA.
- 8) Allow the pellet to dry, then resuspend in 50-75  $\mu\text{l}$  of RNase-Free Water for quantitation. **Do not resuspend the RNA in an EDTA-containing solution** if the RNA will later be enzymatically capped (e.g., with CELLSCRIPT's ScriptCap Cap 1 Capping System).
- 9) Quantitate the RNA by spectrophotometry or fluorimetry. If desired, adjust the concentration of the RNA with RNase-Free Water. The RNA can now be frozen and stored at  $-20^{\circ}\text{C}$  or  $-70^{\circ}\text{C}$ .

II) **Organic Extraction / Ammonium Acetate Precipitation:** Removes all proteins and selectively precipitates RNA, leaving most of the unincorporated NTPs in the supernatant. Note: for this method, the RNA to be purified must be >100 bases in size.

- 1) Adjust reaction volume to 200  $\mu\text{l}$  total using RNase-Free Water (add 179  $\mu\text{l}$  to the reaction).
- 2) Add one volume (200  $\mu\text{l}$ ) of TE-saturated phenol/chloroform. Vortex vigorously for 10 seconds.
- 3) Spin in a microcentrifuge at  $>10,000 \times g$  for 5 minutes to separate the phases.
- 4) Remove the aqueous (upper) phase with a pipette and transfer to a clean tube.
- 5) Add one volume (200  $\mu\text{l}$ ) of 5 M ammonium acetate, mix well then incubate for 15 minutes on ice.
- 6) Pellet the RNA by centrifugation at  $>10,000 \times g$  for 15 minutes at  $4^{\circ}\text{C}$ .
- 7) Remove the supernatant with a pipette and gently rinse the pellet with 70% ethanol.
- 8) Remove the 70% ethanol with a pipette without disturbing the RNA pellet.
- 9) Allow the pellet to dry, then resuspend in 50-75  $\mu\text{l}$  of RNase-Free Water for quantitation. **Do not resuspend the RNA in an EDTA-containing solution** if the RNA will later be enzymatically capped (e.g., with CELLSCRIPT's ScriptCap Cap 1 Capping System).
- 10) Quantitate the RNA by spectrophotometry or fluorimetry. If desired, adjust the concentration of the RNA with RNase-Free Water. The RNA can now be frozen and stored at  $-20^{\circ}\text{C}$  or  $-70^{\circ}\text{C}$ .

Continued on next page.

III) **Organic Extraction / Chromatography / Ethanol Precipitation**: Removes all proteins and unincorporated NTPs from the RNA.

- 1) Adjust reaction volume to 200  $\mu$ l total using RNase-Free Water (add 179  $\mu$ l to the reaction).
- 2) Add one volume (200  $\mu$ l) of TE-saturated phenol/chloroform. Vortex vigorously for 10 seconds.
- 3) Spin in a microcentrifuge at  $>10,000 \times g$  for 5 minutes to separate the phases.
- 4) Remove the aqueous (upper) phase with a pipette and transfer to a clean tube.
- 5) Remove unincorporated NTPs by spin column chromatography.<sup>5</sup> For commercially-available columns, follow the manufacturer's instructions for this step. Recover the RNA in  $\sim 100 \mu$ l.
- 6) Add one-tenth volume (10  $\mu$ l) of 3 M sodium acetate and 2.5 volumes (250  $\mu$ l) of 95% ethanol to the tube, mix well.
- 7) Incubate for 15 minutes on ice.
- 8) Pellet the RNA by centrifugation at  $>10,000 \times g$  for 15 minutes at 4°C.
- 9) Remove the supernatant with a pipette and gently rinse the pellet with 70% ethanol.
- 10) Remove the 70% ethanol with a pipette without disturbing the RNA pellet.
- 11) Allow the pellet to dry, then resuspend in 50-75  $\mu$ l of RNase-Free Water for quantitation. **Do not resuspend the RNA in an EDTA-containing solution** if the RNA will later be enzymatically capped (e.g., with CELLSCRIPT's ScriptCap Cap 1 Capping System).
- 12) Quantitate the RNA by spectrophotometry or fluorimetry. If desired, adjust the concentration of the RNA with RNase-Free Water. The RNA can now be frozen and stored at  $-20^{\circ}\text{C}$  or  $-70^{\circ}\text{C}$ .

IV) **RNA-Binding Purification Column**: Several options are available commercially from multiple vendors. Follow the manufacture's recommended protocol.

- 1) Follow the manufacture's recommended protocol.
- 2) The final resuspension of RNA should be in RNase-Free Water for quantitation. **Do not resuspend the RNA in an EDTA-containing solution** if the RNA will later be enzymatically capped (e.g., with CELLSCRIPT's ScriptCap Cap 1 Capping System).
- 3) Quantitate the RNA by spectrophotometry or fluorimetry. If desired, adjust the concentration of the RNA with RNase-Free Water. The RNA can now be frozen and stored at  $-20^{\circ}\text{C}$  or  $-70^{\circ}\text{C}$ .

## TROUBLESHOOTING

Symptom	Solution
<b>Low yields or less than full-length transcripts</b>	Cleanup the templates to remove any RNase or other contaminants (see Technical Appendix for procedure).
	Verify that ScriptGuard RNase Inhibitor was added to the reaction.
	Extend the incubation time. Do not extend the reaction time beyond 4 hours.
	Increase the template concentration.
	Increase the reaction temperature to 42°C.
	Check for stable secondary structure ( $T_m > 37^\circ\text{C}$ ) in the DNA template which can cause T7 RNAP pausing or dissociation from the template.
<b>Assembled reaction formed an insoluble precipitate</b>	Repeat assembly of the reaction at $>22^\circ\text{C}$ .
<b>White precipitate in reaction buffer</b>	Incubate the reaction buffer at $37^\circ\text{C}$ for 5 minutes then mix thoroughly to dissolve the precipitate.
	Do not store the kit at $-70^\circ\text{C}$ .

## RELATED PRODUCTS

- A-Plus™ Poly(A) Polymerase Tailing Kit
- Cap-Clip™ Acid Pyrophosphatase
- INCOGNITO™ T7-FlashScribe™ Ψ-RNA Transcription Kit
- INCOGNITO™ T7-FlashScribe™ N1meΨ-RNA Transcription Kit
- INCOGNITO™ T7 mScript™ Ψ-mRNA Production System
- INCOGNITO™ T7 mScript™ N1meΨ-mRNA Production System
- INCOGNITO™ SP6 Ψ-RNA Transcription Kit
- INCOGNITO™ T7 5mC- & Ψ-RNA Transcription Kit
- INCOGNITO™ T7 ARCA 5mC- & Ψ-RNA Transcription Kit
- MessageMAX™ T7 ARCA-Capped Message Transcription Kit V2
- ScriptCap™ 2'-O-Methyltransferase Kit
- ScriptCap™ Cap 1 Capping System
- ScriptCap™ m<sup>7</sup>G Capping System
- ScriptGuard™ RNase Inhibitor
- SP6-Scribe™ Standard RNA IVT Kit
- T7-FlashScribe™ Transcription Kit V2
- T7 mScript™ Standard mRNA Production System V2
- T7-Scribe™ Standard RNA IVT Kit

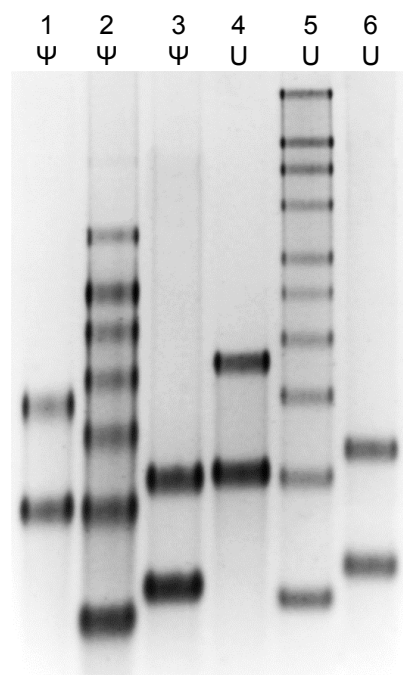
## REFERENCES

1. Karikó, K. et al., (2008) Molecular Therapy 16, 1833.
2. Anderson, B.R. et al., (2010) Nucl. Acids Res. 17, 5884.
3. Karikó, K. et al., (2005) Immunity 23, 165.
4. Karikó, K. and Weissman, D. (2007) Curr. Opin. Drug Discov. Devel. 10, 523.
5. Sambrook, J. et al., (1989) Molecular Cloning: A Laboratory Manual (2nd ed.), New York, Cold Spring Harbor Laboratory Press.

## TECHNICAL APPENDIX

**A. Electrophoretic Mobility of Ψ-RNA in Denaturing Agarose**

Ψ-RNA displays altered mobility during electrophoresis as compared to the comparable U-RNA of identical sequence. This is relevant when comparing Ψ-RNA transcript size to a U-RNA molecular weight ladder.

2% denaturing agarose gel.

Lane 1) Ψ-RNAs, 1,875 & 1,077 bases.

Lane 2) Ψ-RNA Molecular Weight Markers  
(0.5, 1.0, 1.5, 2.0, 2.5, 3.0 & 4.0 kilobases).

Lane 3) Ψ-RNAs, 1,203 & 673 bases.

Lane 4) U-RNAs, 1,875 & 1,077 bases.

Lane 5) U-RNA Molecular Weight Markers  
(0.5, 1.0, 1.5, 2.0, 2.5, 3.0, 4.0, 5.0, 6.0 & 9.0 kilobases).

Lane 6) U-RNAs, 1,203 & 673 bases.

**B. Clean-up of Problematic Templates**

Templates that give low yields or less than full-length transcripts may contain RNase or other contaminants. Such templates usually give better results after the following treatment:<sup>5</sup> See Purification of the Transcription Product Section III, skip step 5.

- 1) Add Proteinase K to 100-200 μg/ml and SDS to 0.5%.
- 2) Incubate for 30-60 minutes at 37°C.
- 3) Extract with an equal volume of a 1:1 mixture of TE-saturated phenol/chloroform.
- 4) Ethanol precipitate.
- 5) Gently remove the supernatant and rinse the pellet with 70% ethanol.
- 6) Resuspend in RNase-Free TE Buffer.



**Research Use Label License**

**INCOGNITO™ kits or other products or services comprising N1-methyl-pseudouridine, pseudouridine, 5-methyl-cytidine and/or other modified nucleotides for use in making RNA or mRNA** that are covered by technologies and IP disclosed in U.S. Patent Nos. 8,278,036 and 9,750,824 and PCT Publ. Nos. WO 2007/024708 and WO/2011/071931 and U.S. and international divisional, continuation, or other patents and patent applications derived therefrom that are owned by the University of Pennsylvania and licensed to CELLSCRIPT ("Products") are sold for research use only.

By purchasing these Products from CELLSCRIPT or an authorized distributor of CELLSCRIPT, Purchaser receives a limited non-exclusive, non-transferable right to use the Products purchased from CELLSCRIPT solely for its own internal laboratory research use (the "Product Use"). Such Product Use expressly excludes any Commercial Use comprising [A] clinical or other use in humans, [B] veterinary, livestock, agricultural or commercial use in animals, and/or [C] manufacture, distribution, importation, exportation, or sale of any products and/or services made using Products, it being understood that [A]–[C] are separately and/or collectively defined as a "Commercial Use" herein. By purchase of Products, Purchaser agrees not to make, import, use, transfer, offer for sale, or sell Products, components of Products, or derivatives of Products, including modified nucleoside-containing RNA or cells made by use of Products, or to provide services, information, or data obtained by the use thereof in exchange for money or other consideration. CELLSCRIPT provides no warranties (statutory or implied) concerning non-infringement of intellectual property rights of any other parties, and all such warranties are expressly disclaimed. Please contact CELLSCRIPT to inquire about a license for a Commercial Use.

The license remains in effect as long as the product is being used in accordance with the license terms. The license is terminated if the user breaches any of the license terms, such as using the product for commercial purposes without proper authorization.

The performance of this product is guaranteed for one year from the date of purchase.

Triton is a registered trademark of Rohm & Haas, Philadelphia, Pennsylvania.

A-Plus, Cap-Clip, CELLSCRIPT, INCOGNITO, MessageMAX, mScript, ScriptCap, ScriptGuard, SP6-Scribe, T7-FlashScribe and T7-Scribe are trademarks of CELLSCRIPT, Madison, Wisconsin.

The Purchaser of this product agrees to the TERMS AND CONDITIONS posted on CELLSCRIPT's website: <http://www.cellscript.com>.

© 2024 CELLSCRIPT, All rights reserved.